


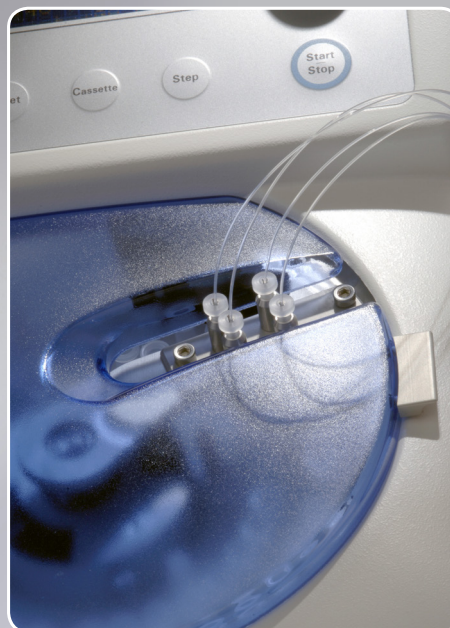


- 温度低至+6°C
- 收集样品量1~1200 L
- 可放置64个收集瓶
- 可同步四通道收集样品

CMA 470 低温微量收集器

The CMA 470低温微量收集器是专为微升体积（ μL ）的微透析样品所设计。此产品电子式的冷却系统可至6°C。也可在密封的收集瓶进行收集，防止样品的蒸发和化学物质的降解，是两个重要的设计方面的考量。可收集每管样品至瓶底小到1 μL 的样品。可放入64个300 μL 或40个2 mL收集瓶。CMA 470有不同的收集固定支架，可同时进行四通道的分步收集。

低温控制容器内有可旋转的收集瓶卡槽。收集固定支架是垂直悬空的，不会影响收集的流速，CMA 470收集器不需要连接电脑就能正常运行，但CMA 470还是配有一个RS232和一个USB端口可与其他仪器或电脑进行连接或进行升级。

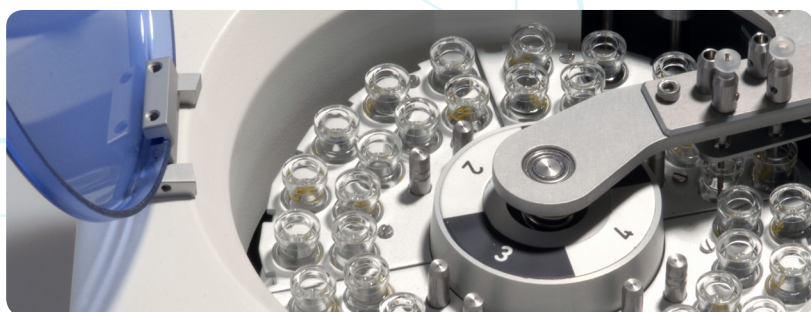



技术指标

收集量:	1 μ L - 1.2 mL
收集样品数:	64 x 300 μ L, 40 x 2 mL
密封盖:	一次性使用
低温:	可至 +6 °C
温度精准度:	$\pm$ 1.5 %
冷却能力:	室温降至15 °C以下
收集模式:	时间, 分, 秒
工作环境温度:	10 – 32 °C
显示:	LCD 4 x 20 characters, 背光模式
存储:	可存储设置
外部连接:	通用数字输入/输出
计算机连接:	RS 232串行接口和USB
电压:	100 – 240 VAC, 50 – 60 Hz
功率:	外部适配器最大55 W
尺寸:	222 x 279 x 142 (167) mm (W x D x H)
重量:	约3.8 kg

设计用途: 设计用于研究和工业应用, CMA 470冷冻部分收集器未批准用于临床。

订货信息

	Ref.No
CMA 470 Refrigerated Fraction Collector	CMA8002770
including:	
4 Cassettes for Small Vials	
Holder for Single Cannula	
Holder for Dual Cannulae	
Holder for Quadruple Cannulae	
4 Cannulae for Tubing	
4 Cannulae for Septa	

配件

	Ref.No.
Vials, Plastic 300 μ L (1000/pkg)	CMA7431100
Caps Plastic re-sealing 1000/pkg*	CMA7431102
Vials, Glass 300 μ L (500/pkg)	CMA7431007
Caps/Seals Non-Reclosing, Small, 1000/pkg*	CMA7432175
Crimper, Small	CMA7432017
Cassette, Small Vials, Plastic	CMA8320010
Cassette, Large Vials, Aluminum	CMA8320008
Cannula for Tubing	CMA8002999
Cannula for Septa	CMA003000

* 不能在收集时使用